


The Entellus™ PathAssist™
Light Fiber™ is a versatile
illumination tool that offers
real-time confirmation
during endoscopic
sinus procedures

An economical confirmation tool designed to work with the
XprESS™ Multi-Sinus Dilation Tool and Entellus Sinus Cannula



The PathAssist Light Fiber is a versatile tool that provides economical illumination and confirmation during endoscopic sinus procedures

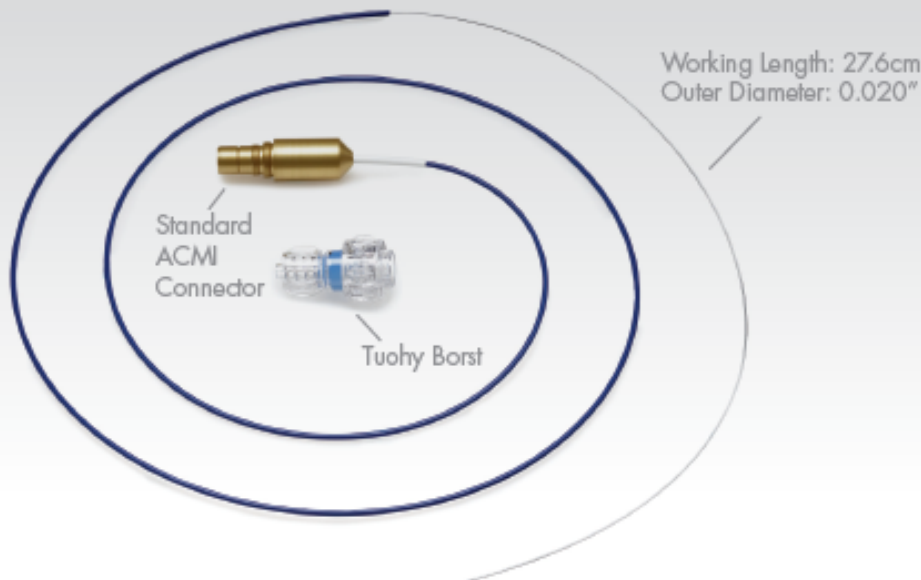

- Provides transillumination to confirm correct entry to the frontal recess and location of the frontal sinus
- Delivers additional directable light to illuminate dark sinus pathways
- Useful during FESS, nasal endoscopy, and balloon dilation procedures
- Small diameter device suitable for office use on awake patients
- Alternative or additive confirmation device to CT - image guidance systems or fluoroscopy
- 4x more light output than other transillumination devices
- For single use only

Compatible with multiple devices

The Entellus Light Fiber is designed for use with the provided Entellus Sinus Cannula, the XprESS™ Multi-Sinus Dilation Tool, malleable suctions, and other compatible devices.

For added convenience, the Light Fiber has an ACMI connection and can easily be connected to Wolf or Storz light cables with adaptors.

Light Fiber


Compatible with working lumen instruments with an internal lumen diameter $\geq 0.035"$ and a length $\leq 27.5\text{cm}$



High-intensity output for transcutaneous confirmation

Balloon inflation port

An ideal tool for real-time illuminated confirmation with the XprESS Multi-Sinus Dilation Tool and other compatible devices

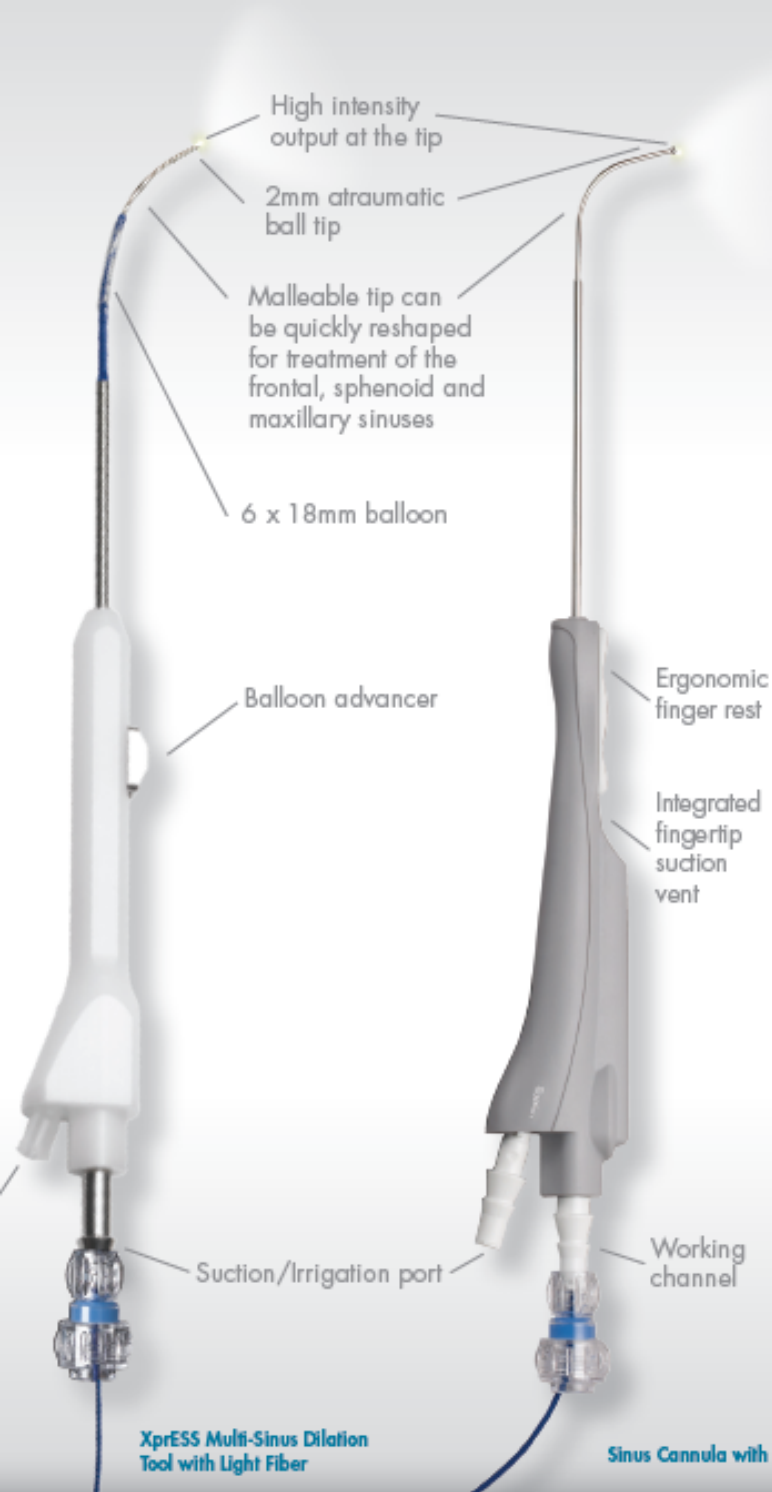

Designed for use with the XprESS Multi-Sinus Dilation Tool

When used with the Entellus XprESS Multi-Sinus Dilation Tool, the Light Fiber enables balloon inflation with transcutaneous light confirmation of tool position. After dilation, simply deflate the balloon and remove the Light Fiber from the XprESS device to resume suction or irrigation.

Economical, disposable illumination confirmation with the Entellus Sinus Cannula

Use the shapeable Entellus Sinus Cannula with Light Fiber to gently confirm entry into a sinus prior to FESS, suction, irrigation, or dilation procedures. Also, useful as a screening tool to determine accessibility of a sinus prior to treatment. After confirming correct location in the sinus, simply remove the Light Fiber and proceed with irrigation or suction.

- Easy-to-use, intuitive design provides tactile feedback
- Well-suited for sinus procedures in the office or OR
- Compatible with multiple light source connections



The Sinus Cannula can be used for suction and irrigation during sinus procedures

The Entellus PathAssist Light Fiber is a simple, practical solution for in-office illumination

The Entellus Light Fiber is a versatile addition to the Entellus PathAssist family. An economical option for physicians who want transcutaneous illumination to confirm location of the sinus, the Light Fiber is designed for efficient use with the provided Sinus Cannula, the XprESS Multi-Sinus Dilation Tool and other compatible devices.

To learn more about the Entellus Light Fiber, XprESS Multi-Sinus Dilation Tool and the advantages of in-office balloon sinus dilation procedures, visit www.EntellusMedical.com

About Entellus Medical™

Entellus Medical delivers innovative, high-quality, and minimally invasive therapeutic solutions to healthcare providers and their patients who suffer from sinusitis.



To order: Please call 866-620-7615

CATALOG NO.	DESCRIPTION
LF-100	PathAssist Light Fiber (includes Sinus Cannula)
JD-100	XprESS Multi-Sinus Dilation Tool

XprESS Indications for Uses:

To access and treat the frontal recesses, sphenoid sinus ostia, and maxillary sinus ostia/ethmoid infundibula in adults using a trans-nasal approach. The bony sinus outflow tracts are remodeled by balloon displacement of adjacent bone and paranasal sinus structures.

Light Fiber Indications for Uses:

To locate, illuminate within, and transilluminate across nasal and sinus structures in adults aged 18 and over.

Please see products' Instructions For Use for a complete listing of contraindications, warnings, precautions, and adverse events.

Caution: Federal (USA) law restricts this device to sale by or on the order of a physician.

©2011 Entellus Medical. PathAssist, Light Fiber and XprESS are trademarks of Entellus Medical, Inc. 2125008 rA 11/11



Entellus Medical
6705 Wedgwood Court North
Maple Grove, MN 55311
866-620-7615
www.EntellusMedical.com